



Non-Latex Surgical Gloves

A SECOND SKIN... THAT DOESN'T IRRITATE THE FIRST



Delivering Unsurpassed Performance and Protection for Every Hand

The GAMMEX® Non-Latex product line offers the most complete range of non-latex options, taking allergy prevention beyond latex with advanced protection against chemical allergic reactions for both patients and healthcare workers.

Providing Comprehensive Protection for Healthcare Workers and Patients

Ansell uses proprietary technologies and formulations to create a wide range of non-latex options that are easy to use and provide protection and safety for any allergy profile and surgical need.

- Polyisoprene options offer **latex-like comfort** to healthcare workers suffering from latex allergies and type IV allergies caused by diphenylguanidine (DPG)
 - **DPG-free polyisoprene gloves** include:
 - GAMMEX® PF IsoDerm®
 - GAMMEX® PF IsoDerm® Sensitive
 - GAMMEX® PF IsoDerm® Underglove
- An **accelerator-free neoprene** option, GAMMEX® PF DermaPrene®, offers extra protection to healthcare workers who have type IV allergies to chemical accelerators

The Performance You Expect From GAMMEX®

Like all GAMMEX® gloves, the Non-Latex product line combines unmatched comfort and sensitivity with strength. Each glove is produced on Ansell's proprietary formers for a natural fit that reduces hand fatigue and stress.


GAMMEX® Non-Latex gloves deliver confirmed protection by:

- Meeting or exceeding international standards
- Delivering strong resistance to viral permeation, verified by extensive bacteriophage test results
- Offering a protective barrier for handling mixed bone cement product*
- Providing high resistance to alcohol-based hand disinfectants

The line includes options with SureFit Technology™, which eliminates cuff roll-down and prevents fluid from entering the glove.

* Provide a secure barrier only when standard preparation and use directions for bone cement are followed.

GAMMEX® Non-Latex range minimizes the risks of type IV allergies and skin irritations, but some chemicals contained in GAMMEX® Non-Latex gloves may cause allergic reactions or irritations in some users.

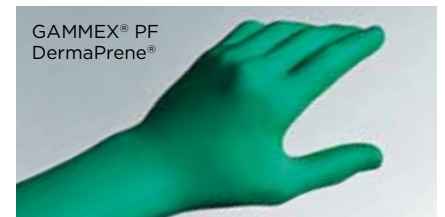



Non-Latex Surgical Gloves

ANSELL – DELIVERING ON OUR COMMITMENT TO MINIMIZE ALLERGIES WITHIN THE HEALTHCARE ENVIRONMENT

A Full Line of Non-Latex Products Offering Best-in-Class Hand Protection

The broad GAMMEX® Non-Latex product line meets an array of surgical needs while reducing risk and providing protection for patients and healthcare workers with latex and chemical allergies.

Special Features

- A cetylpyridinium chloride-free polymer coating with hydrophilic and hydrophobic components offers easy donning and changing with wet and dry hands
- SureFit Technology™, available on GAMMEX® polyisoprene gloves, eliminates cuff roll-down and provides surgeons, nurses, and other healthcare professionals with superior protection during operating room procedures
- An anti-slip smooth finish facilitates better instrument grip in both damp and dry environments
- Ansell's proprietary anatomic formers with curved fingers provide comfort with a natural fit

Quality Standards

- Extended product testing for enhanced quality assurance
- Statistically tested for pinholes to G1; AQL of 1.0 after packing
- Conforms to EN 455, AS/NZ 4179, and ASTM D3577
- Manufactured in compliance with ISO 9001:2008, ISO 13485:2003, and ISO 14001:2004

Glove Material and Formulation

	COLOR	MATERIAL	ACCELERATOR FREE	DPG-FREE	POWDER-FREE	FREE OF ANIMAL DERIVATIVES
GAMMEX® PF DermaPrene®	Green	Polychloroprene; no natural rubber	Yes	Yes	Yes	Yes
GAMMEX® PF IsoDerm®	White	Polyisoprene; no natural rubber		Yes	Yes	Yes
GAMMEX® PF IsoDerm® Sensitive	White	Polyisoprene; no natural rubber		Yes	Yes	Yes
GAMMEX® PF IsoDerm® Underglove	Green	Polyisoprene; no natural rubber		Yes	Yes	Yes

Dimensional Specifications

	LENGTH (typical average values, mm)	
GAMMEX® PF DermaPrene®	≥305	(Size 5.5-Size 9)
GAMMEX® PF IsoDerm®	≥300	(Size 5.5-Size 9)
GAMMEX® PF IsoDerm® Sensitive	≥300	(Size 5.5-Size 9)
GAMMEX® PF IsoDerm® Underglove	≥300	(Size 5.5-Size 9)

Thickness (Single)

	FINGER	PALM	CUFF
(typical average values, mm)			
GAMMEX® PF DermaPrene®	0.185	0.175	0.16
GAMMEX® PF IsoDerm®	0.23	0.215	0.175
GAMMEX® PF IsoDerm® Sensitive	0.19	0.17	0.14
GAMMEX® PF IsoDerm® Underglove	0.19	0.17	0.14

Ordering Information

	ORDER CODES SIZE 5.5 TO 9	PACKAGING
GAMMEX® PF DermaPrene®	353201 to 353208	25 pairs per dispenser box; 4 boxes per case
GAMMEX® PF IsoDerm®	353611 to 353618	40 pairs per dispenser box; 6 boxes per case
GAMMEX® PF IsoDerm® Sensitive	353601 to 353608	
GAMMEX® PF IsoDerm® Underglove	353621 to 353628	

Neither this document nor any other statement made herein by or on behalf of Ansell should be construed as a warranty of merchantability or that any Ansell product is fit for a particular purpose. Ansell assumes no responsibility for the suitability or adequacy of an end user's selection of gloves for a specific application.



Ansell, ® and ™ are trademarks owned by Ansell Limited or one of its affiliates. © 2011 Ansell Limited. All rights reserved.

WWW.ANSELL.EU

Medical Solutions

